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Protect the environment! Do not throw electrical items in the bin. Take them to the nearest recycling centre so that they can be destroyed or recycled.

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UTILISATION

NOTICE

1. Please use the voltage of the machine identified.
2. After using and before cleaning must be unplugged.
3. Not in moist places.
4. Children or people with disabilities, should be paid special attention.
5. Pregnant women, medical patients do not use.
6. If there is a problem, do not repair by yourself, do not continue to use it, you should contact your dealer as soon as possible.

EMERGENCY TREATMENT

1. If this product is abnormal, please turn off the power switch, then unplug the power plug.
2. If internal mechanical failure, please send to our company service station around, our company provides a dedicated service.
3. If the power line is damaged, shall replace it by the manufacturer or service or similar qualified personnel so that avoid the danger.

FAULT DIAGNOSIS

Fault	Reason	Solution
Motor not working	Plug in not well	Plug in the power & plug well
	Switch is not in place	Check the switch

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
ACCESSORIES INTRODUCTION

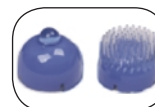Massage head

Three different shapes of the massage cover design, when massaging it can concentrate strength and easy to operate, free to move and massage in all parts of the body.

Three-combinations

Three different combinations mode, the user can replace them at will, let you in the enjoyment of happiness and happiness of enjoying,
The locking strap design make the product can firmly control in the hands and arbitrary all parts of the body.

Massage strength adjustment

Continuously Variable slide type push button switch, from comfortable relaxation to excited refreshing, makes you feel warm in use.

One key switch design

The switch has three gears, the first gear to start the vibration, the second gear to start the vibration with red light, the third gear to turn off the internal power supply. It makes the operation more convenient and the use more intimate.

HUMAN BODY ACUPOINT DIAGRAMS

The body from neck to waist, there are about 650 acupoints, and every acupoints is closely related to certain organs or physical condition. By massage, kneading and other techniques on the body, to promote blood circulation, improve metabolism and other effects. Massage hammer make you get the most comfortable pleasure and health of fun in a relaxed way!

SPECIFICATION AND DIAMETERS

Product name:	body massager
Item No.:	LY-618
Rated Power:	25 W
Frequency:	50 Hz
Rated Voltage:	AC 220V
Working Times:	15 mins

The company will continuously improve the various models of the performance and appearance. If there is a change without notice ,please understand.

USER'S MANUAL ATSU ELECTRIC MASSAGER Ref. 15193220

Thank you for choosing our product. We hope you enjoy the device.

Before using the device for the first time, please read the instructions for use very carefully and keep these safe along with the guarantee certificate, receipt and if possible the box and the inner packaging. If you give the device to third parties, you should also give them the instructions for use.

FUNCTION AND FEATURES


1. Split structure design, three combinations, light and handy, convenient operation , let you enjoy the pleasure of massage anytime ,anywhere.
2. The main part of the massager can be used alone, relying on the locking strap fixed in the palm to make the operation more simple, more humane.
3. Massager's handle is convertible. When the handle is installed, you can massage the common parts of the body; When the handle is anti installed, you can massage the special parts of the body; so as to achieve the effect of body massage.
4. Continuously Variable speed regulation can effectively stimulate the nerves and muscles, to reach comfortable and relaxed massage effect.
5. Point massage helps relieve muscle pain, let you feel relaxed and full of vigor.

OPERATION METHOD

1. Plug in the power plug.
2. Press the power switch once to start the massage function. Self-liking,regulating The speed control button to choose the appropriate massage strength.
3. Press the power switch twice to start the massage function with red light.
4. To stop the massager, just press the power switch again.
5. After using , Please remember to unplug.
6. Every massage is appropriate for 15 minutes, if massage over 30 minutes, excessive massage

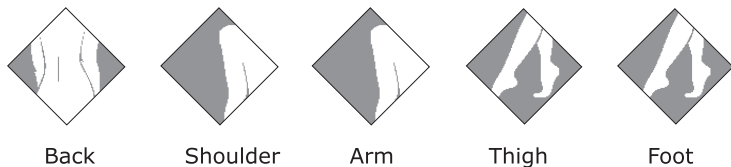May over-stimulate the muscles and nerves, and thus be counterproductive.

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