NATURE &
DECOUVERTES

REVEIL-LAMPE DE LUMINOTHERAPIE 2-EN-1

2-IN-1 LIGHT THERAPY WAKE-UP LAMP

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

WARNING

- Do not use this appliance in wet surroundings (e.g. in the bathroom or near a shower). Do not immerse appliance in water or any other liquid.
- Before use, ensure that all packaging materials are removed and that there is no visible damage to the appliance or accessories. When in doubt, do not use the appliance and contact your dealer or the customer service address provided.
- The adaptor contains a transformer. Do not cut off the adaptor to replace it with another plug, as this can cause a hazardous situation.
- This appliance is only intended for use on the human body.
- Close supervision is necessary when this appliance is used by, or near children or invalids.
- Children should be supervised to ensure that they do not play with the appliance.
- There is no need for any special training for this appliance.
- Keep the cord away from hot objects and open flames.
- Connect the appliance only to the mains voltage listed on the type plate.
- Always remove the mains plug and let the appliance cool off before you touch it.
- Do not touch the appliance with wet hands while it is plugged in; do not allow any water to be sprayed onto it. The appliance must be operated only when it is completely dry.
- Protect the appliance from strong impacts.
- Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- If the adaptor, cord or the appliance is damaged, contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.
- Disconnection from the power supply network is only guaranteed when the mains part is unplugged.
- This appliance is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- This appliance requires no calibration and preventive inspection.
- This appliance is not repairable and contains no user serviceable parts.
- The operator shall not touch the appliance and the patient simultaneously.

- The appliance is protected against unauthorized use. Do not modify this appliance without authorization of the manufacturer.
- If this appliance is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.

GENERAL NOTES

Important

- This appliance is not intended for commercial or clinical use, but only for individual, private household use.
- If you are taking medication such as pain relievers, antihypertensives or antidepressants, consult your doctor before using this appliance.
- People with retinal diseases as well as diabetics need to be examined by an ophthalmologist before using this appliance.
- Please do not use if you have an eye disease such as cataracts, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humour.
- If you have health concerns of any kind, contact your general practitioner!

EMC INFORMATION

Warning

- Warning: Use of this appliance adjacent to or stacked with other device should be avoided because it could result in improper operation. If such use is necessary, this appliance and the other device should be observed to verify that they are operating normally. If the appliance does not light on or flickering, it means that the appliance is not operating normally.
- Warning: Do not change any cable without authorization of the manufacturer. Use of accessories and cables other than those specified or provided by the manufacturer of this appliance could result in increased electromagnetic emissions or decreased electromagnetic immunity of this appliance and result in improper operation.
- Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.
- Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm to any part of this appliance (GLO3), including power cable. Otherwise, degradation of the performance of this equipment could result.
- If the appliance does not light on or flickering due to the electromagnetic disturbances, do not use it and try to place the appliance in other location.
- For maintaining the basic safety and essential performance with regard to electromagnetic disturbances for the expected service life, please precaution to take if the use location is near (e.g. less than 1,5 km from) from AM, FM or TV broadcast antennas.

MAINTENANCE AND REPAIRS

Important

- Avoid exposing the screen to direct sunlight, spotlight or any other heat source.
- Do not use alcohol or ammonia base liquid to clean the screen. The cleaning information provided in the section "Cleaning".
- Please do not attempt to repair the appliance yourself. This could result in serious injuries. Failure to observe this regulation shall void the warranty. If this appliance is damaged in any way or it is not working, it should not be operated. Please try to fix the problem according to FAQ table. If the problem is still not resolved, please contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.

INSTRUCTIONS FOR USE

Congratulations on the purchase of your new therapy lamp. The product will offer many years of reliable usage if handled correctly. Please keep your receipt, these instructions and the original packaging safe in case of future requirement.

Please read the instructions carefully to gain the most benefit from your new therapy lamp, although the usage guidelines in this document do not replace the advice of your qualified medical advisor.

INTENDED USE

Medical intended use of this appliance is to treat the symptoms of SAD (seasonal affective disorder) or winter depression. This therapy lamp is intended to boost people's energy levels, to adjust circadian rhythms and make people feel more awake. This therapy lamp is safe and medically proven to provide relief from the symptoms of SAD.

HOW DOES THE THERAPY LAMP WORK?

In autumn and winter, the seasons with the least sunlight because of shorter days, the lack of sunlight and a problem with certain brain chemicals stops the hypothalamus working properly. The lack of light is thought to affect the production of the hormone melatonin (that makes you sleepy).

This therapy lamp uses bright lights to simulate the sunlight during the darker autumn and winter, whilst eliminating the dangerous UV rays contained in sunlight. The additional light helps your brain to reduce to production of hormone melatonin, adjust the levels of melatonin that can help the symptoms of SAD or winter depression.

SEASONAL AFFECTIVE DISORDER

Seasonal Affective Disorder (SAD) as known as winter depression or winter blues. It may cause people:

- Feel low, listless and fatigued
- Need for more sleep
- Have trouble concentrating
- Less active than normal
- Lack energy

Symptoms occur during the same season each year. Most people with SAD start to have symptoms during winter and feel better in summer.

CONTRAINDICATION

The following conditions may be contraindication for use of this appliance:

- suffer(ed) from severe depression
- recent eye surgery or a diagnosed eye condition whereby doctor has advised you to avoid bright light on specific medication (i.e. certain antidepressants, psychotropic drugs or malaria tablets);
- have high blood pressure with vitreous hemorrhage;
- skin sensitive to light, such as systemic lupus erythematosus

Always consult your doctor before your start using this appliance if you have any of the above conditions.

Possible side effects:

- Headaches
- Eyestrain
- Sleeping problems
- Nausea

Measures to avoid, minimize or alleviate these side effects:

- Increase the treatment distance
- Avoid using the appliance during the late evening or before sleep
- Stop using the appliance for a few days to make the side effects disappear and then try again.

SAFETY AND PERFORMANCE STANDARDS

This therapy lamp has been designed and manufactured to conform to the highest standards of safety and performance including electromagnetic compatibility (EMC). The therapy lamp conforms to the applicable requirements of the following:

- This appliance meets the requirement of Medical Device Directive 93/42/EEC for medical appliances.
- The appliance is complied with the following emission and immunity standard, and test level:

Emissions Test Standard	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/flicker emissions IEC 61000-3-3	Comply

GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY							
	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
Radiated RF IEC61000-4-3 (Test specification for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	385	380 - 390	TETRA 400	Pulse modulation 18 Hz	1,8	0.3	27
	450	430 - 470	GMRS 460, FRS 460	FM \pm 5 kHz deviation 1 kHz sine	2	0.3	28
	710	704 - 787	Bande LTE 13, 17	Pulse modulation 217 Hz	0,2	0.3	9
	745						
	780						
	810	800 - 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, Bande LTE 5	Pulse modulation 18 Hz	0,2	0.3	28
	870						
	930						
	1720	1700 - 1990	GSM 1800; CDMA 1900; GSM 1900; DECT; Bande LTE 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	0,2	0.3	28
	1845						
	1970						
	2450	2400 - 2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, Bande LTE 7	Pulse modulation 217 Hz	0,2	0.3	28
	5240	5100 - 5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	2	0.3	9
	5240						
5785							


Note: UT is the a.c. mains voltage prior to application of the test level.

SIGNS AND SYMBOLS

The following symbols appear in this manual, on the therapy lamp or on its accessories:

	Warning instruction indicating a risk of injury or damage to health		Read operating instructions before use
	Safety note indicating possible damage to the appliance/ accessory		Caution: Consult accompanying documents
	Specifies serial number of the therapy lamp		Class II Equipment
	Daylight ON/OFF		Date of Manufacture
	Authorized Representative in the European Community		Manufacturer
IP21	Ingress Protection Rating		Keep dry

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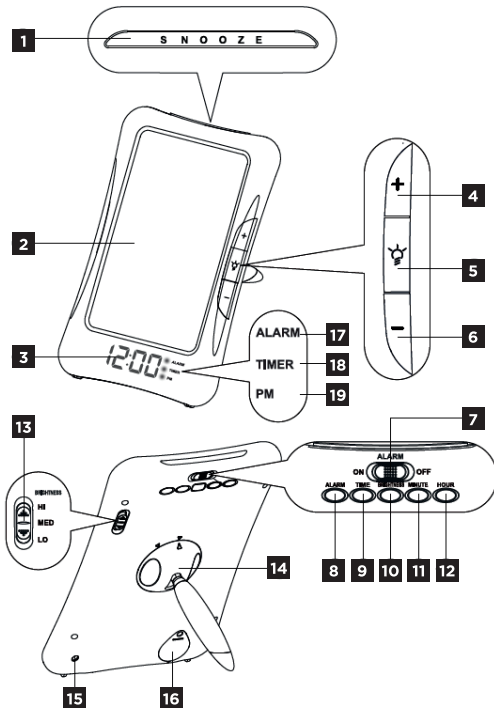- 1 Snooze
- 2 Screen
- 3 Display panel
- 4 +
- 5 Light
- 6 -
- 7 Alarm ON/OFF
- 8 Alarm Set
- 9 Time Set
- 10 Clock Brightness
- 11 Minute
- 12 Hour
- 13 Wake Up Brightness
- 14 Removable stand
- 15 DC input
- 16 Clock backup battery
- 17 Alarm indicator
- 18 Timer indicator
- 19 PM indicator

INSTALLATION

1. Insert the end of stand into the hole on the back of appliance.
 - a. Insert the stand into the cavity on the back of the unit, with the triangular mark on the circular part of the stand aligning with the 'UNLOCK' triangular mark on the back housing.
 - b. After that, turn the stand in clockwise direction until the mark on the stand meets the 'LOCK' triangular mark.
 - c. There should be a 'click' sound and then you will feel snug.
 - d. Loosen the handle by rotating it in anti-clockwise direction.
 - e. Move the handle vertically to reach your desired angle.
 - f. Tighten the handle by rotating it in clockwise direction.
2. Place the appliance on a solid surface approximately 20 cm from your eyes, the light must bathe your face from the side.
3. Connect the appliance only to the mains voltage listed on the type plate..

Note

- Push the plug completely into the mains socket.
- Ensure that there is a socket near the place where the appliance will stand.
- Lay the mains cable in such a way that no one can trip over it.
- Do not rotate the handle too tight to prevent the handle mechanism damage.

HOW TO OPERATE

Before use, please check the appliance for damage or faults. If you notice any damage or faults on the appliance, do not use it and contact customer services or your supplier.

Daylight

1. Press "Light" button at the face side of the device to turn on the daylight therapy lamp (the lamp).
2. At the same time, red timer will flash on the display panel.
3. Select the therapy time (15, 30, 45, 60 minutes or no selection) by pressing the "Light" button repeatedly.
4. The countdown of therapy time will start once the red timer stopped flashing. The remaining time will be indicated.
5. The time can be reset to 15, 30, 45, 60 minutes or no selection by pressing the "Light" button during the countdown. If you don't select the therapy time, the appliance will turn OFF after 4 hours.
6. There are 5 levels of brightness which can be selected by pressing the "+" or "-" button. (The default level is 3)
7. The lamp will automatically turn on when countdown finished (shown as "00").
8. Or, press "Light" button and held for 2 seconds to turn off the lamp. When the lamp is turned off, the red timer will be turned off as well.
9. You can read, write, make telephone calls or work while you use the appliance. When using the appliance, try to maintain the recommended distance of 17-27 cm between your eyes and the appliance. During the darkest period of the year, repeat the treatment for at least 7 consecutive days, or longer, depending on your individual requirements.

**Usage guidelines are as follows:
(in brightness level 5)**

Distance	Duration
17 cm	Approx. 1 hours
27 cm	Approx. 2 hours

Note

- Do not look directly into the light continuously! Some eye pain or mild headaches may be experienced during initial sessions but these should pass. If in doubt, just lower your treatment times until you feel comfortable again.
- Do not leave the appliance unattended in the duration.

SET-UP THE CLOCK

1. Press the “Time Set” button at the back of the device and hold for 3 seconds.
2. After that, the clock will flash on the display panel at the face side of the device.
3. Adjust the clock by pressing the “Hour” and “Minute” buttons repeatedly. For example, to set the clock as 7:00am, press “Hour” button to shown “7” and press “Minute” button to show “00”. Blue light of “PM” will be shown when setting time in p.m.
4. Confirm the clock by pressing the “Time Set” button again.
5. The clock setting will also be confirmed automatically if no operation is performed within 10 seconds and the digits will stop flashing.

ADJUST THE CLOCK BRIGHTNESS

1. Adjust the brightness of display panel by pressing the “Clock Brightness” button (at the back of the device) repeatedly.
2. There are 4 levels of brightness including High, Middle, Low and Invisible. The default brightness is High.

SUNRISE LIGHT

Set-up the Alarm

1. Slide the “Alarm ON/OFF” button (at the back of the device) to “I” and thus to turn ON the alarm.
2. Press and hold the “Alarm Set” button for 3 seconds.
3. After that, the alarm time will flash on the display panel at the face side of the device.
4. Adjust the alarm time by pressing the “Hour” and “Minute” buttons repeatedly. For example, to set an alarm at 7:00am, press “Hour” button to shown “7” and press “Minute” button to show “00”. Blue light of “PM” will be shown when setting time in p.m.
5. Confirm the alarm time by pressing the “Alarm Set” button again.
6. The alarm setting will also be confirmed automatically if no operation is performed within 10 seconds and the digits will stop flashing.

USING THE SUNRISE LIGHT

1. Slide the “Wake Up Brightness” button at the back of device to adjust the brightness levels of the sunrise light.
2. There are 3 levels of brightness (from Hight to Low) for selection. For beginner, would suggest to start from Low level (the mildest sunrise light effect).
3. The yellow LED sunrise effect will start at 30 minutes before the alarm time set. For example, if the alarm clock is set at 7:00am, the sunrise light will start at 6:30am. The alarm will ring (a beep sound) at 7:00am.
4. Press the “Snooze” button at the top of the device to stop the alarm.
5. Or, just leave the alarm rings for 4 minutes. It will stop automatically. However, this temporary stop of alarm will not disturb the sunrise light.
6. After 9 minutes, the alarm will ring again.
7. Slide the “Alarm ON/OFF” button (at the back of the device) to “O” and thus to turn OFF the both alarm and sunrise light.

Note: The sunrise light must not be considered for medical purpose.

STORAGE

When this appliance is not being used, it should be stored in a dry location out of reach of children. Do not keep the appliance connected to a power source. Follow the storage instructions provided in the section “Technical Specification”.

CLEANING

Warning

- Before you clean the appliance, always switch it off, unplug it and let it cool off.
- Do not touch the appliance with wet hands.
- Clean using a damp cloth, to which you can apply a little detergent if necessary.
- Do not let water run into the appliance or spill water onto the appliance.
- Cleaning can be done once per month.

REPLACING CLOCK BACKUP BATTERY

- This appliance is equipped with a battery backup circuit that will maintain the time and alarm time settings in the event of a temporary power interruption. The battery backup circuit requires one CR2032 battery, included.
- When replacing the backup batteries be sure to have the unit connected to the AC outlet to maintain all of your settings, otherwise it will be necessary to reset the time and alarm time after the battery is replaced. Follow steps 1 - 3 below to install the replacement batteries.
 1. Remove the battery cover located on the back of the appliance and remove the exhausted CR2032 battery.
 2. Insert 1 fresh CR2032 battery in the backup battery compartment with the “+/-” ends facing as shown.
 3. Replace the battery compartment cover.

Note:

- CR2032 battery is included in the unit. Please pull out the battery pull tab before you use the appliance.
- Make sure you dispose of battery properly (see your local regulations).
- Do not charge or recharge the battery.
- Do not install backward.
- Do not dispose of in fire or heat above 75°C (167°F).
- Do not expose contents to water

FAQ

Problem	Possible causes	Solutions
Clock does not light on	The brightness of clock selected in invisible	Press the clock brightness button to adjust the brightness
	No Power	Connect the main plug correctly
	No Power	The main cable is defective. Contact customer services or your retailer
Daylight does not light up	Light button switched off	Switch on the light button
	LEDs have reached the end of its service life or faulty.	If you need to have the appliance repaired, contact customer service or an authorized dealer.
Alarm has not sound	Speaker faulty	If you need to have the appliance repaired, contact customer service or an authorized dealer.
Sunrise does not light up	Setting error	Please recheck the alarm setting according to the instruction
	LEDs have reached the end of its service life or faulty.	If you need to have the appliance repaired, contact customer service or an authorized dealer.

TECHNICAL SPECIFICATION

- Model: GLO3
 - Dimension (WHD): 260 x 170 x 23mm
 - Weight: 600g
 - Lighting elements: LEDs
 - Power: DC12V, 1A
 - Light intensity: Around 10,000 lux in 10 cm
 - Mains connection: 12V DC, 1.0A, only in adaptor model: HP24A-120100-AdV or HP24A-120100-AdB
 - Operating Environment: 0 to +40 oC, and 30% to 90% relative humidity
Atmospheric pressure: 700-1060hPa
 - Storage and Transport Environment: -10 to +60 oC, and 30% to 90% relative humidity
Atmospheric pressure: 700-1060hPa
 - Pression atmosphérique : 700hpa-1060hpa
 - Product Classification:
 1. Class II equipment
 2. No applied part
 3. Protection against ingress of water: IP21
 4. Not category AP/APG equipment
 5. Mode of operation: Continuous operation.
 - Maximum length of power cord: 2000 mm
 - Location of emission aperture: Screen
 - Treatment area: Around 1800 cm²
 - Risk Group of light source equipment
Group 1 in blue light classified by IEC 60601-2-57:2011
 - Product Service Life: 5 years
- Subject to technical modification.

Maximum output of optical radiation for the “GLO3”:

Radiance Output	Risk group classified by IEC 60601-2-57:2011	Max. value
EUVA: Eye UV-A	Exempt Group	6,69×10 ⁻⁵ W·m ⁻²
ES: Actinic UV skin & eye	Exempt Group	1,9×10 ⁻⁵ W·m ⁻²
EIR: Infrared radiation hazard exposure limits for the eye	Exempt Group	2,34×10 ⁻² W·m ⁻²
EH: Thermal hazard for the skin	Exempt Group	1,34×10 ¹ W·m ⁻²
LB: Blue light	Risk Group 1	3,86×10 ² W·m ⁻² ·sr ⁻¹
LR: Retinal thermal	Exempt Group	8,16×10 ¹ W·m ⁻² ·sr ⁻¹

Important

If this appliance is damaged in any way or it is not working and you cannot find the answer to your question in FAQ. Please contact customer service or the retailer for assistance.

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