
**NATURE &  
DECOUVERTES**

**SYSTEME PAR PRESSION D'AIR**  
AIR PRESSURE SYSTEM  
SISTEMA DE PRESIÓN DE AIRE  
SYSTEEM MET LUCHTDRIUK  
SISTEMA POR PRESSÃO DE AR  
**Réf. 15225530**

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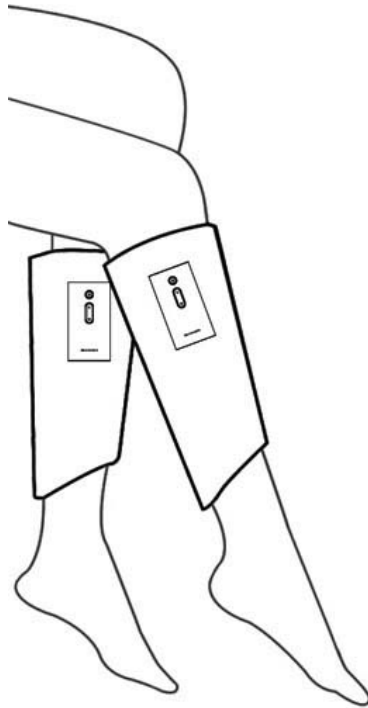# PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

## PRODUCT INTRODUCTION

Our leg massager is simple while fashionable in design easy to operate 'massaging the leg through the way of air bag extrusion, can effectively alleviate edema/muscle ache/soreness after exercise has good effect on leg shaping. It will become your ideal family massage master.

This product mainly through the multi-control air bag repeatedly, orderly inflation and deflation, imitating the real masseur's soft deep tissue pressing technology. Through hot compress to improve local blood circulation and nutrition supply, enhance muscle vitality, keep calf muscles in a relaxed state.



Easy to operate Three working modes, three gear vibration force adjustment and three heating temperatures. It can be used alone or on both legs.

- Product Name: Leg Massager Product model: D01-CM006
- Rated power: 5W
- Rated voltage: 5V
- Battery capacity: 2500mAh
- Default timing: 15 mins/time
- Charging time: 4hours, running times: 10 times-15times



## INSTRUCTIONS

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- This appliance contains batteries that are non-replaceable.
- **WARNING:** For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. (ZL-005L0501000EU02)
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

1. When you use it for the first time, if it is low-battery, please charge the device before using. While low-battery, the power indicator will flicker
2. Pull off the magic sticker of device, allowing the product to be fully unfolded.
3. Please bond the massage instrument to the lower leg, ensure that the power on key facing the right, and adjust to appropriate tightness. (Note: the product can be used while clothing on)
4. Long press  power key to turn on the device, short press the power key to switch mode, by default mode one, a total of three modes in cycle.

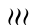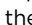5. Working state, short press  Key to turn on the hot compress model, a total of three sets of hot compresses (laboratory measurement data, please note the low temperature scald).
6. Working state, short press the  key to turn on vibration mode, a total of three strength for adjustment.
7. The device automatically shuts down for 15 minutes of working time.
8. While USB charging, the power indicator is incrementally flicking, four lights all on when battery fully charged. The product does not work under charging condition, while charging, there is no effect to press any key.

## KEY DESCRIPTION
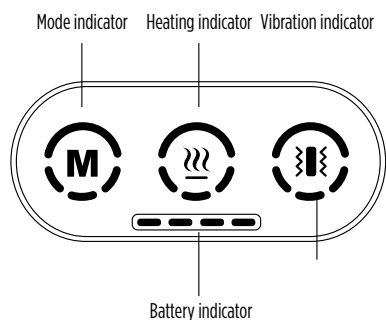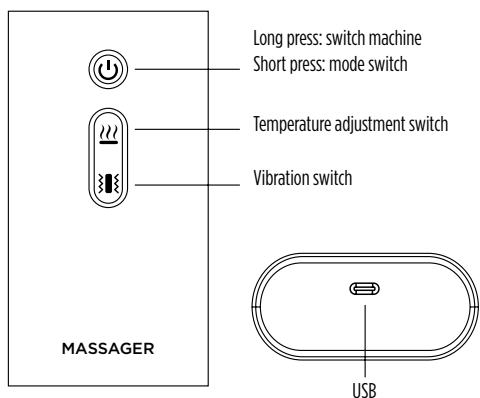
## CLEANING MAINTENANCE

- The main engine should be well protected and cleaned. Wipe with a dry and soft cloth. Avoid moist cloth or a cleaning agent. Avoid water leakage.
- Product with cloth area does not supportive for water washing, if needed. Please wipe with alcohol pad or clean wet towel. (Gasoline or gasoline diluent are forbidden)
- Please keep this product out of the reach of children, avoid moist environment, high temperature, direct sunlight and water. If not using the device for a long period. Keep the battery charged every 2months.
- Before abandoning the appliance power off the appliance and take out the lithium battery, the battery must be disposed according to local regulations.

## AVERTISSEMENTS

- Do not use while it is charging
- While using. Turn off the power first if intending to stop using.
- Do not use while driving in inflammable explosive or damp environment.
- Please avoid sharp instruments contacting with this product, so as not to damage the product.
- This product contains lithium battery, before abandoning the appliance, power off the appliance and take out the lithium battery the battery must be disposed according to local regulations, the battery can not be treated as domestic waste, do not throw the battery into fire.
- Do not use during sleep.
- Only professional staff allowed to disassemble or repair this product.
- Do not use when power cord is damaged.
- Do not use when water of other liquid leakage.
- Do not use when the product is abnormal.
- Stop using immediately if you have any discomfort.
- Do not place this product.

## TROUBLESHOOTING

Problem	Possible causes	Solutions
The device working too Low strength or can not feel it at all	Device not started?	Veillez le mettre en marche
	Improper tightness?	Veillez le repositionner, l'ajuster.
	Low battery?	Veillez charger l'appareil.
The device can not be started	Is the battery full ?	Chargez complètement l'appareil avant de l'utiliser.
There is no heating in the massage machine	The heating function is not turned on?	Long press M key to turn on

## IMPORTANT

IT SHOULD BE USED WITH CAUTION OR UNDER THE GUIDANCE OF MEDICAL DOCTORS IN CASE OF THE FOLLOWING SITUATIONS.

- Pregnant or postpartum women
- Suffering from bone disease
- Skin allergy (e.g. scald or burn)
- Patients with severe heart disease and abnormal blood pressure
- Patients with frost epilepsy
- The body has implanted devices (e.g. pacemakers)
- People with osteoporosis

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