**NATURE &
DECOUVERTES**

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

PRECAUTIONS

Danger. Please never use the product with the following medical electronic devices:

1. Heart pacemakers and other implantable medical electronic devices
2. Artificial heart and lungs and other life-supporting medical electronic devices;
3. Electrocardiographs and other medical electronic devices;

* These medical electronic devices may cause malfunctions and can be life-threatening.

Warning:

- This appliance contains batteries that are only replaceable by skilled persons.
This appliance contains batteries that are non-replaceable.
- Use only with ZL-005L0501000EU02 supply unit
- Warning: For the purpose of recharging the battery, only use the detachable supply unit provided with this appliance.
- Regarding the instruction for information regarding charging refer to the section "CHARGING THE DEVICE".
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- If the patient applies high-frequency surgical device and the product on the body at the same time, the electrode of the product may cause burns and damage to the product.
- If the product is used near (1M) a short-wave or micro-wave therapy device, the output of the product may be unstable.
- The electrode should prevent contact with wounds and scars during treatment.
- When you use the product, don't place the electrode near the chest heart, otherwise the risk of heart fibrillation will increase.
- The product is recommended for use by one user. When it is used by different users, clean the electrode on the product before using it.
- When you use the product, stick the electrode tightly and evenly to the skin.
- The electrode of this product has a heat-generating function. When the hot compress function is applied on the same part for a long time (more than 15 minutes), it may cause low temperature burns.
- If you feel unwell during or after use, stop using it and seek medical attention.
- In the following cases, consult your doctor before use, otherwise it may cause accidents or cause physical discomfort.

1. Patients under medical treatment and feel discomfort.
2. Patients with abnormal heart and brain nerves
3. Patients whose body temperature exceeds 38 degrees (fever period).
4. Patients with infectious diseases.
5. Patients with skin perceptual disorders or skin abnormalities.
6. Patients under me dictation.
7. Patients with abnormal blood pressure.
8. Patients who are not free to express their wishes.
9. Patients with sensory impairment due to high peripheral circulatory disorders such as diabetes.
10. Patients with skin allergies to stainless steel electrodes.

• Usage under the following conditions is prohibited:

1. Malignant tumor patients
2. Patients with bleeding tendency
3. Pregnant women, lying-in women
4. Acute suppurative inflammation patients

5. Heart pacemaker users
6. Patients with metal implant
7. Near the heart

- Don't apply it near the heart, above the neck, head, mouth, genitals, or skin disorders and don't treat the soles of the feet simultaneously. In addition, don't compress the internal organs when using. Otherwise, an accident or physical discomfort may occur.
- Don't use with other medical devices or more than 2 products. Otherwise, nausea or physical discomfort may occur.
- Don't use for purposes other than the ones described in the instruction manual. Otherwise, accident or malfunction may occur.
- Don't use the electrode with a chipped or damaged surface. Otherwise, accidents or malfunctions may occur.
- Don't disassemble, repair or modify this product without authorization. Otherwise, malfunction or electric shock may occur.

Note:

- If you feel unwell due to an abnormality in this product, stop using it immediately.
- If your body or skin feel strange, consult your doctor and follow the doctor's advice.
- If you want to move the pad or place it other on other person durine massage, be sure to turn off the power before moving it. Otherwise you may receive a strong electric shock.
- Don't treat while having electronic devices. Otherwise, the electronic device may malfunction(for example, the electronic watch disorder).
- If abnormal operation occurs, immediately switch off the power supply and remove it. Otherwise, high temperatures or short circuits may occur and lead to an accident. At the same time contact customer service for after-sales advice.
- Don't plug or unplug the USB charging cable with wet hands. Doing so may result in electric shock or accident.
- Don't use it for children or people who have no expressiveness, otherwise it may cause an accident or cause discomfort.
- Don't use it in a place with high humidity, such as a bathroom, or use it while bathing. Otherwise, you may receive a strong electric shock and damaged the product.
- Don't use it when you are sleeping. Otherwise, it may cause the main body to malfunction. In addition, the pad may move to other parts and cause discomfort.
- Don't use while driving a car, otherwise it may cause physical discomfort due to excessive stimulation, resulting in an accident.

- Don't stick the electrode to the affected area for a long time otherwise it may cause skin inflammation.
- When massaging, don't use the metal parts of belts, necklaces, etc to touch the electrode, otherwise you may receive a strong electric shock.
- Don't apply the product on the same place for long time (more than 30 minutes), otherwise it may cause muscle fatigue at the massage site to cause physical discomfort.
- Please wear this product correctly and use it on the cervical spine. Wear it symmetrically. Otherwise, accidents or physical discomfort may occur.
- When unplugging the USB charging cable from the body do not pull on the cable, but rather, pull the cable plug and pull it out. otherwise, it may result in electric shock or damage to the product.
- The product electrode has a surface heating function for users who are not sensitive to heat or subject to heat allergies.
- Please consult a doctor before use.
- When the product is working for a long time, the surface will be hot, Which is normal.
- When the remote control battery is exhausted, replace it with a new battery of the same model as soon as possible. Battery model: CR2032
- When replacing the remote control battery, do not install the pole reversibly.
- After the remote control battery is exhausted dispose of it according to the relevant regulations in your place of residence

Suggestion :

- Please make sure that the host button is working properly.
- When you use it for the first time, or if the battery loses power after the product has been left idle for a long time, please perform the charging operation according to the instruction manual to ensure that it can work normally.
- When disposing of this product, please dispose of it in accordance with the relevant environmental protection regulations of the city, otherwise it may cause environmental pollution.
- It is recommended for people aged 16 and above.
- Don't use or store this product in the following locations:
 - Places exposed to direct sunlight
 - Places subject to high temperature and humidity
 - Places with water
 - Dusty places
 - Places near the fireworks
 - Places susceptible to vibration shocks
 - Strong electromagnetic places.

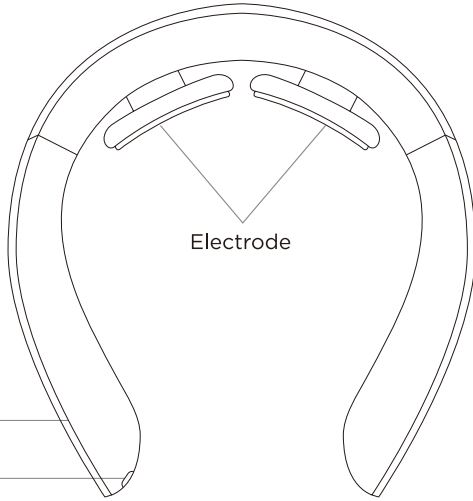
CHARGING THE DEVICE

1. Before first use, charge the product fully for about 2 hours.
2. The digital display flashes regularly, indicating that the device is charging.
3. The digital display stays lit once the battery is fully charged.

Note:



- The battery can be charged at any time and at all power levels.
- This product cannot work while charging, it should be used after being disconnected from the power supply.
- If unused for a prolonged period of time, recharge the battery completely every three months.

PRODUCT DIAGRAM


Button areas

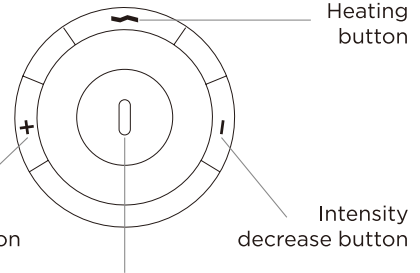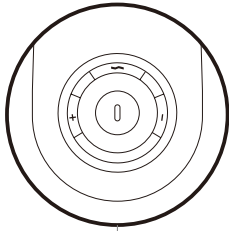Charging interface


Switch Mode button

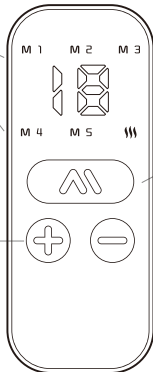Screen area(5 modes display, Force digital display, Heating mode display)

Massage mode and heating mode button

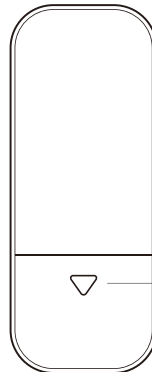Intensity button (increase and decrease)

OPERATION METHOD

The massager can be divided into Ordinary type and Honorable type:

The Ordinary style has only one massager, while The Honorable style also has a remote control and a pair of multi-function electrodes.

1. Operation mode before USE:

- Before use, first wet the skin of the cervical spine with a wet towel, and then put the product around the neck, make the electrode close to the skin of the cervical spine, adjust the position of the product, host and adjust the appropriate massage intensity and mode. then put the product around the neck, make the electrode close to the skin of the cervical spine, adjust the position of the product, host and adjust the appropriate massage intensity and mode.

2. Startup & Shutdown:

- Long press the host « 0 » pbutton for two seconds, turn on the host (with the buzzer “ beep ” sound and blue light), enter the default mode with the default level (0).
- Long press host « 0 » for two seconds when the machine is on, hear two “beep” sounds and turn off the products.

3. Mode choice (currently 5 models)

- After the product host starts up, short press the host « 0 » as the mode key, short press the mode key or the mode key on the remote control to switch the mode.
- There are five modes. After each press, it will send out “ beep ” sound and switch to higher mode. The five modes can be switched in cycle (Note: M1 mode, M2 mode, M3 mode, M4 mode, M5 mode). Each time a mode is switched, the intensity returns to the default 0 level.

4. Gear selection

- “ + ” on the remote control or the product host is the intensity increase button. After each press, it will send out “beep” sound and switch to higher mode. The intensity is divided into 15 levels (0, 1, 2... 15). When the intensity is at level 15, When pressing” + “, it will send out two “beep” sound with same intensity at level 15.
- « - » on the remote control or the product host is the intensity decrease button. After each press, it will send out “beep” sound and switch to lower mode. The intensity is divided into 15 levels (0,1,2 ... 15).when the intensity is at level 0,it will send out two “beep” sound with same intensity at level 15 after pressing « - ».
- There is a display screen on the remote control, which will display the mode selection and gear (intensity) selection as well as the heating function. After pressing any key on the remote control, the massager will respond with a “beep” sound (when the massager is on), and the display will be turned off after 5 seconds.

5. Heater selection

- « ☺ » on the product host or « ☺ » the remote control is the heating button. Press the heating key to turn on the heating function of the host machine with “beep”, and the constant temperature of 40-42 °C will be reached gradually within 30 seconds.
- Press « ☺ » the product host tock for 2 seconds to turn off the heating button or to press « ☺ » the remote control tock for 2 seconds. Press the tock button to turn off the heating button to stop the heating function.

6. Timing function

- The host is scheduled for 15 minutes. When the host is started up, it will automatically shut down after 15 minutes. It is recommended to use it for no more than 30 minutes every day.
- When the product host is on, if the product is not normally worn and used, it will automatically shut down in 1 minute.

7. Multi-functional electrode plates can be used

- The Honorable style is equipped with a set of multi-functional electrodes, which can be inserted by the charging port of the massager to enjoy the same massage experience as the current state of the massager (no heating mode).

8. Low voltage alarm & charging function

- Low voltage alarm: when the host is short of power (battery voltage is 3.5v), the blue light of the indicator flashes quickly (frequency: 3 times per second), you can insert one end of micro-usb cable into the host and the other end into the mobile phone Charger (output is 5V 1A/2A), green light breathing flash when charging, green light always on when full.

CURRENT MAINTENANCE

- When cleaning and maintaining this product, please turn off the product host, clean the electrode with a wet paper towel, and wipe and disinfect with a little medical alcohol (75% alcohol concentration) when it is especially dirty.
- When daily maintenance, please use neutral detergent and soft cloth to wipe, do not use corrosive detergent, do not put into the water to soak or rinse.
- When daily maintenance, please use neutral detergent and soft cloth to wipe, do not use corrosive detergent, do not put into the water to soak or rinse.

TROUBLESHOOTING

- If you have problems using the product, please refer to the following tips and solutions.
- If the product is still in a fault condition or cannot be used, turn off the power.
- Users should not disassemble or repair the product by themselves. Dismantling or repairing by themselves will cause damage to the product and cause the invalidity of warranty. The company will not be responsible for any consequences caused thereby.

Problems	Possible cause/treatment
Failure to start	<ul style="list-style-type: none">• The internal rechargeable battery is running low• The remote control battery is dead.• Please confirm whether the host product switch has been pressed.
Failure to charge or no charge light	<ul style="list-style-type: none">• Please check if the USB cable and host and phone power adapter are connected.
Normal display but no stimulation	<ul style="list-style-type: none">• Does the pad contact the skin directly? Adjust the position so that it is close to the skin.• Is there hair or clothing, jewelry, etc. between the pads and the skin? Move away or remove items.• Is the intensity level at 0? Increase the massage intensity.
Tingling sensation	<ul style="list-style-type: none">• Is the massager's electrode and neck skin closely contacted? adjust the position so that it's stick to the skin.• Is the massage time too long? It is recommended to massage for 15 minutes at a time, twice a day.• Is the massage intensity too high? Adjust to low level .• Is the skin allergic to the electrode? Those with severe allergies need to be desensitized first.• Is the climate cold and the skin dry? Use a wet towel to moisten the skin and use it.
Weak simulation	<ul style="list-style-type: none">• Is the massager's electrode and neck skin closely contacted? adjust the position so that it's stick to the skin and adjust intensity.• Is the electrode dirty or oily? Clean the electrode surface.• Is the battery running low? Please charge if the battery is low.• Is massage intensity level too low? Increase massage intensity.

PRODUCT PARAMETERS

Main material: PC+silicone

Battery (product): Lithium battery 650mAh.

Battery(controller): Lithium cell CR2032

Input Voltage: DC5V/1A

Power: 5W

The static current: <100uA

The company reserves the right to improve the product and the manual, if there is any change without notice, the appearance of the product in kind shall prevail. This book has been carefully checked, but there may be printing mistakes or misunderstanding of the content, please understand, the company reserves the right to explain.

Precautions during handling of batteries

When used correctly, lithium batteries provide a safe and dependable source of power. However, if they are misused or abused, leakage, venting or in extreme cases, explosion and/or fire can result.

- a) Keep batteries out of the reach of children

 WARNING	
KEEP OUT OF REACH OF CHILDREN Swallowing can lead to chemical burns, perforation of soft tissue, and death. Severe burns can occur within 2 hours of ingestion. Seek medical attention immediately.	

- b) Do not allow children to replace batteries without adult supervision.
- c) Always insert batteries correctly with regard to polarity (+ and –) marked on the battery and the equipment.
- d) Do not short-circuit batteries.
- e) Do not charge batteries.
- f) Do not force discharge batteries.
- g) Do not mix new and used batteries or batteries of different types or brands.
- h) Exhausted batteries should be immediately removed from equipment and properly disposed of.
- i) Do not heat batteries.
- j) Do not weld or solder directly to batteries.
- k) Do not dismantle batteries.
- l) Do not deform batteries.
- m) Do not dispose of batteries in fire.
- n) A lithium battery with a damaged container should not be exposed to water.
- o) Do not encapsulate and/or modify batteries.
- p) Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.
- q) Remove batteries from equipment if it is not to be used for an extended period of time unless it is for emergency purposes.

Disposal of Batteries

Batteries may be disposed of via communal refuse arrangements provided no local rules to the contrary exist.

During transport, storage and handling for disposal, the following precautions should be considered:

- a) Do not dismantle batteries.
- b) Do not dispose of batteries in fire except under conditions of approved and controlled incineration.
- c) Store collected batteries in a clean and dry environment out of direct sunlight and away from extreme heat.
- d) Store collected batteries in a well-ventilated area.
- e) Do not mix collected batteries with other materials.
- f) Protect battery terminals.

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- 1) This appliance contains batteries that are only replaceable by skilled persons.
Instructions for appliances containing non-replaceable batteries shall state the substance of the following: This appliance contains batteries that are non-replaceable.
- 2) Use only with ZL-005L0501000EU02 supply unit
- 3) Warning: For the purpose of recharging the battery, only use the detachable supply unit provided with this appliance.
- 4) Specifications of USB charger:
Ratings: 5.0V --- 1.0A

Information published	Value and precision	Unit
Manufacturer's name or trade mark, and address.	Shenzhen Zhongli Power Technology Co., Ltd. Bldg. F 2-4 Layer, Li Hao Yongda High-Tech Park Gongming ST, Guangming New District, 518105 Shenzhen, Guangdong, P.R China	
Model number	ZL-005L0501000EU02	
Input voltage	100-240	V
Input AC frequency	50/60	Hz
Output voltage	DC5	V
Output current	1	A
Output power	5	W
Average active efficiency	≥75,88	%
Efficiency at low load(25%)	≥ 75,15	%
No-load power consumption	≤0,1	W

5) Battery Warnings:

- The supply terminals are not to be short-circuited.
- Attention should be drawn to the environmental aspects of battery disposal. Don't throw used batteries in dust bin. Please contact your retailer in order to protect the environment.
- The batteries (batteries installed) shall not be exposed to excessive heat such as sunshine, fire or the like.
- If battery leakage occurs, avoid contact with skin. Isolate the leaking battery in a sealed plastic bag and put it in the scrap following the instructions environmental protection.
- If contact with the skin, mucous membranes or eyes, rinse thoroughly with water and contact your doctor or ophthalmologist.
- WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

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